



Connecting the Docking Station

To connect the docking station, follow these steps:

- 1 To select the output on the docking station, do one of the following:
 - To connect the docking station to a computer, skip this step.
 - To connect the docking station to a video monitor, use the video monitor switch on the bottom of the docking station.

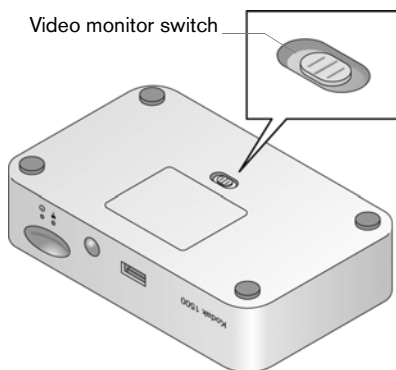


- 2 To connect the docking station, do one of the following:
 - To connect to a computer, connect the computer USB cable to the docking station.
 - To connect to a video monitor, select the correct video cable, VGA, S-video, or AV-out, and connect the cable to the docking station.


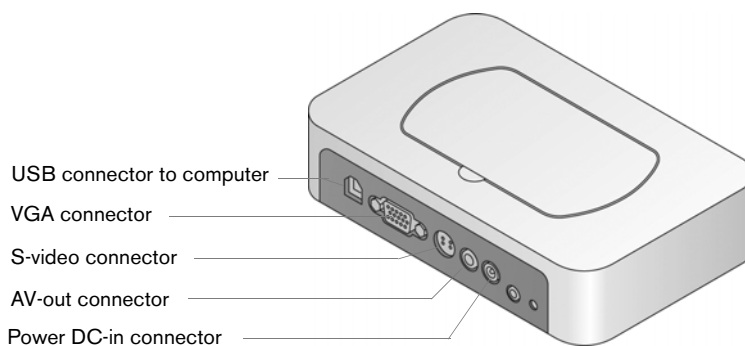


3 To store your images, do the following:

- To connect to a computer, skip this step.
- To connect to a video monitor, insert the SD card into the docking station.

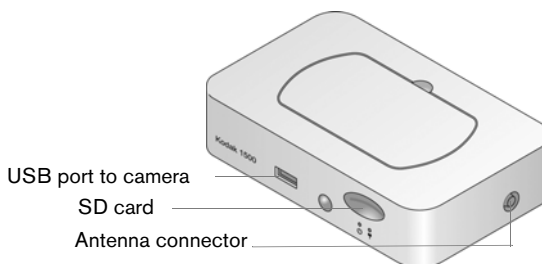

4 Connect the power source to the docking station.

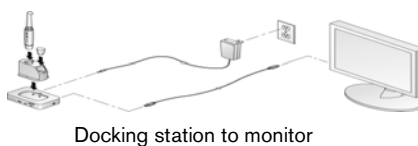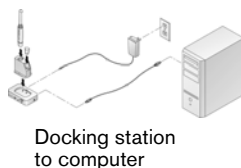

5 If you are using a wireless camera, follow these steps:

- Screw the antenna into the docking station.
- Remove the top cover of the docking station and connect the charge station.

